

Wound dressing combined with secondly dressing material (Ex. gauze)

Plus moist TOP

Plus moist TOP is a new type of wound dressing developed on new product concept of how to use. Plus moist TOP provides a non-adherent function to various conventional dressings by their self-regulating absorption of the exudate fluid. Plus moist TOP is ideal for the wide variety of wounds with various amount exudates, because Plus moist TOP is able to take effective adjustment of absorption quantity by selecting the secondly dressing material which Plus moist TOP combines with.

Features:

- 1) Non-adherent (It does not dry wound surface).
- 2) Able to take effective adjustment of absorption quantity by selecting the secondly dressing material which is integrated with.
- 3) Unique Self-Regulating Absorption of wound exudate.
- 4) Superior breathability appropriate for dermatological therapy.
- 5) No rashes or maceration (skin deterioration).
- 6) Able to cut for any size and shape.



Direction for use:



Table. Absorption quantity of various secondly dressing materials

Secondly dressing material	Absorption quantity [g/100cm ²]
Direct to clothes	—
Gauze (8-ply)	13
Sheet for pets	25
Incontinence pad	40
Diaper	70

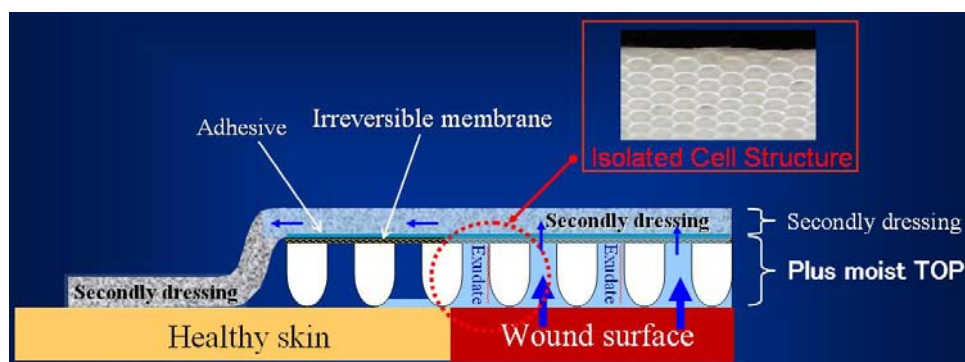


Fig. Function as a self-regulating absorption of the exudate

Product Line

Code	Product Name	Size	Package	Remark
TOP	Plus moist TOP	250 × 200mm	20pcs / bag	

WARNING

Do not reuse. Single use only. Do not use on infected wound which indicates signs such as fever, pain, redness, itching, rash or burning etc. Diabetics and those with poor circulation or wound healing problems should not use this product.